

8353 Genu Direxa wraparound

8356 Genu Direxa

DE Gebrauchsanweisung	5
EN Instructions for use	8
FR Instructions d'utilisation	12
IT Istruzioni per l'uso	15
ES Instrucciones de uso	19
NL Gebruiksaanwijzing	23
SV Bruksanvisning	26
DA Brugsanvisning	30
NO Bruksanvisning	33
FI Käyttöohje	37
PL Instrukcja użytkowania	40
HU Használati utasítás	44
CS Návod k použití	47
HR Upute za uporabu	51
BG Инструкция за употреба	54
TR Kullanma talimatı	58
ET Kasutusjuhis	61
RU Руководство по применению	65
ZH 使用说明书	69

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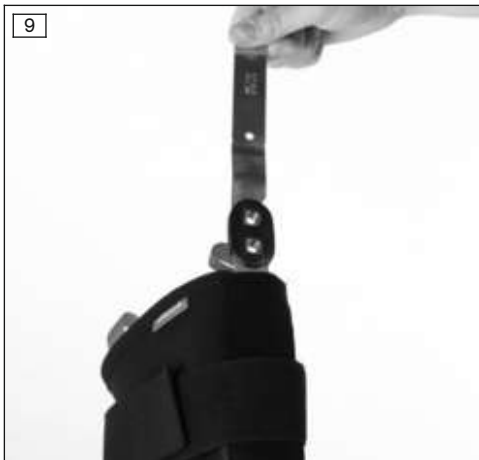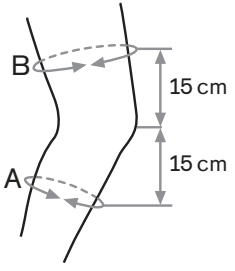INFORMATION

Date of last update: 2021-01-13

- ▶ Please read this document carefully before using the product and observe the safety notices.
- ▶ Instruct the user in the safe use of the product.
- ▶ Please contact the manufacturer if you have questions about the product or in case of problems.
- ▶ Report each serious incident related to the product to the manufacturer and to the relevant authority in your country. This is particularly important when there is a decline in the health state.
- ▶ Please keep this document for your records.

These instructions for use provide important information on the fitting and application of the 8356 Genu Direxa and 8353 Genu Direxa wraparound knee joint braces.

2 Intended use

2.1 Indications for use

The brace is intended **exclusively** for orthotic fittings of the lower limbs and **exclusively** for contact with intact skin.

The brace must be used in accordance with the indications.

2.2 Indications

- Knee pain with mild to moderate ligament laxity and/or mild to moderate feeling of instability
- Posttraumatic and postoperative instability and/or feeling of instability
- Gonarthrosis with initial to moderate degenerative ligament instability and/or slight to moderate feeling of instability (mainly involving the cruciate ligament)
- Rheumatoid arthritis
- Patellofemoral pain syndrome (e.g. chondromalacia patellae, chondropathia patellae, retropatellar osteoarthritis)

The indication must be determined by the physician.

2.3 Contraindications

2.3.1 Absolute Contraindications

None known.

2.3.2 Relative Contraindications

The following indications require consultation with a physician: skin diseases/injuries, inflammation, prominent scars that are swollen, reddening and hyperthermia of the fitted limb/body area; lymphatic flow disorders, including unclear soft tissue swelling distant to the body area to which the medical device will be applied; sensory and circulatory disorders in the leg area.

2.4 Mechanism of Action

The brace improves neuromuscular function, proprioception and circulation. It supports sensorimotor functions and circulation, promotes the resorption of edemas and can relieve pain.

3 Safety

3.1 Explanation of warning symbols

 CAUTION	Warning regarding possible risks of accident or injury.
--	---

 NOTICE	Warning regarding possible technical damage.
--	--

3.2 General safety instructions

 CAUTION

Contact with heat, embers or fire

Risk of injury (such as burns) and risk of product damage

- ▶ Keep the product away from open flames, embers and other sources of heat.

 CAUTION

Reuse on other persons and improper cleaning

Skin irritation, formation of eczema or infections due to contamination with germs

- ▶ The product may be used by one person only.
- ▶ Clean the product regularly.

NOTICE**Improper use and changes**

Change in or loss of functionality as well as damage to the product

- ▶ Use this product with care and only for its intended purpose.
- ▶ Do not make any improper changes to the product.

NOTICE**Contact with oils, salves, lotions or other products that contain oils or acids**

Insufficient stabilization due to loss of material functionality

- ▶ Do not expose the product to oils, salves, lotions or other products that contain oils or acids.

4 Handling

INFORMATION

- ▶ The daily duration of use and period of application are generally determined by the physician.
- ▶ The initial fitting and application of the product must be carried out by qualified personnel.
- ▶ Instruct the patient in the handling and care of the product.
- ▶ Instruct the patient to see a physician immediately if any exceptional changes are noted (e.g. worsening of the complaint).

4.1 Selecting the Size

- 1) Measure the circumference of the thigh 15 cm above the center of the patella.
- 2) Measure the circumference of the lower leg 15 cm below the center of the patella.
- 3) Determine the size of the brace (see size chart).

4.2 Fitting and Application

NOTICE**Use of a worn or damaged product**

Limited effectiveness

- ▶ Before each use, check the product for functional reliability and for possible wear or damage.
- ▶ Do not continue using a product that is no longer functional, or that is worn or damaged.

4.2.1 Adaptation

The joint bars are fitted individually by trained qualified personnel.

4.2.2 Application

⚠ CAUTION**Use with heat-sensitive skin**

Skin irritation due to overheating

- ▶ Do not use the product in case of a known heat allergy.
- ▶ Do not continue to use the product if skin irritation occurs.

⚠ CAUTION**Incorrect or excessively tight application**

Risk of local pressure and constriction of blood vessels and nerves due to improper or excessively tight application

► Ensure that the product is applied properly and fits correctly.

8356 Genu Direxa

- > Open all hook-and-loop closures on the brace.
- > The patient is sitting on the front edge of a chair.
- 1) Hold the brace on both sides with the donning loops and pull it over the knee until the patella cutout is centered on the patella (see fig. 1 and see fig. 2).
- 2) Close the hook-and-loop straps first below and then above the patella (see fig. 2 and see fig. 3).

8353 Genu Direxa wraparound

- > Open all hook-and-loop closures on the brace.
- > The patient is sitting on the front edge of a chair.
- 1) Wrap the brace around the knee from behind (see fig. 4).
- 2) Center the patella cutout over the patella.
- 3) Fasten the hook-and-loop closure below the patella cutout while applying tension. Make sure that the hook-and-loop sections are placed over one another (see fig. 5).
- 4) Fasten the hook-and-loop closure above the patella cutout while applying tension. Make sure that the hook-and-loop sections are placed over one another (see fig. 5).
- 5) Fasten the straps first below and then above the patella. Make sure that the hook-and-loop sections are placed over one another (see fig. 6 and see fig. 7).

4.3 Cleaning

NOTICE

Use of improper cleaning agents

Damage to the product due to use of improper cleaning agents

- Only clean the product with the approved cleaning agents.

- 1) Pull the black strap upwards to remove the bracket over the joint bar (see fig. 8).
- 2) Pull the joint bars out of the brace (see fig. 9).
- 3) Fasten all hook-and-loop closures.
- 4) Recommendation: Use a laundry bag or net.
- 5) Wash the brace in warm water at **40 °C (104 °F)** with standard mild detergent. Do not use fabric softener. Rinse well.
- 6) Allow to air dry. Do not expose to direct heat sources (e.g. sunlight, stove or radiator).
- 7) Reinsert the joint bars.

INFORMATION: The labelling on the joint bars must be visible from the outside.

- 8) Pull the black strap first upwards and then downwards to position the bracket over the joint bar (see fig. 8).

5 Disposal

Dispose of the product in accordance with national regulations.

6 Legal Information

All legal conditions are subject to the respective national law of the country of use and may vary accordingly.

6.1 Local Legal Information

Legal information that applies **exclusively** to specific countries is written in the official language of the respective country of use after this chapter.

6.2 Liability

The manufacturer shall be liable in the event that the product is used in accordance with the descriptions and instructions in this document. The manufacturer will not assume liability for damage caused by disregarding the information in this document, particularly due to improper use or unauthorized modification of the product.

6.3 CE Conformity

The product meets the requirements of Regulation (EU) 2017/745 on medical devices. The CE declaration of conformity can be downloaded from the manufacturer's website.

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